



reevaa **H**

Garvanic ion Scalp Massager

사용설명서 / USER MANUAL

- Please read the manual before using it.  
Please keep it where you can see it at any time after reading it.
- This product is intended for beauty, so use or release for medical treatment is prohibited by law.
- As this product is based on the premise that healthy skin is used, the user is responsible for all the troubles when using it in the area being treated by a doctor.

# PRECAUTIONS ON USE

## Precautions for safety

1. With the battery cover (magnetic ground) held by hand, the brush must touch the scalp to produce galvanic ions.
2. If there is any water left after use, wipe off the water completely and store it after drying.
3. Insert the battery in the right direction when replacing it.
4. After using the machine, remove any foreign substances from the brush and keep it clean.
5. Do not shock or throw the device.
6. Children and the elderly should not use it by themselves.

## Don't use this kind of person.

1. A person suffering from a skin hemorrhage or skin disease or an allergy.
2. Surgery, injury or inflammation in the area of use.
3. Patients with hypertension, hemophilia and heart disease.
4. Pregnant women and menstrual menstruation.
5. A person who uses in-body medical devices, life-support electronic devices such as artificial cardiopulmonary systems, and equipped medical electronic devices such as electrocardiographs.


## Product features

- 1 Iontophoresis by using Galvanic ion
- 2 Long-term use
- 3 Precision Design for User Safety
- 4 Vibration Massage Function
- 5 The brush(consumable item) can be replaced for hygienic use

## The name of each part





## Components

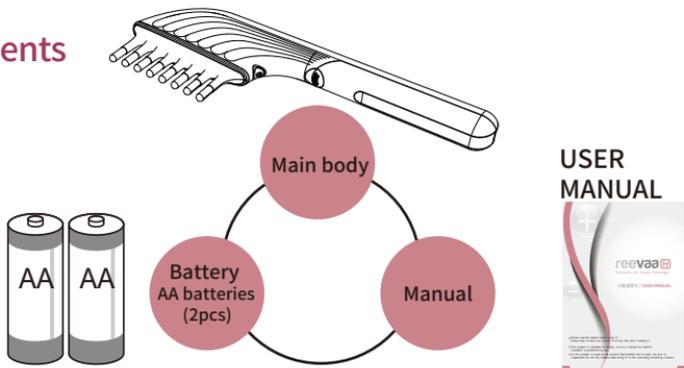

# HOW TO USE

## How to Use

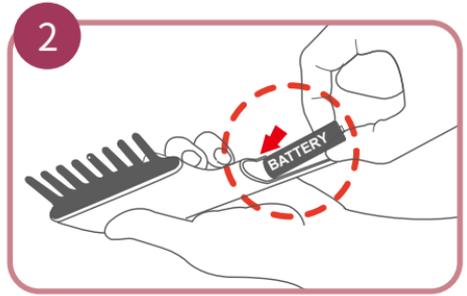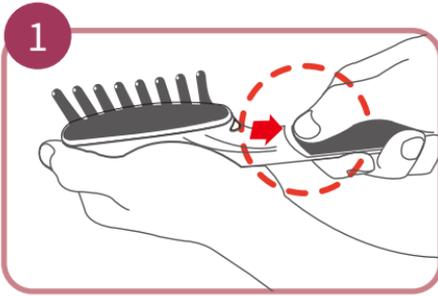

1. Open the battery cover and
2. Put the 2 AA batteries in the direction as shown.  
Insert and close the battery cover.
3. Apply scalp/hair dedicated products, to the scalp/hair.
4. Press the mode button to set it to  $\oplus$  ions or  $\ominus$  ions for the product.

\* If there is no separate mark such as product, set to  $\ominus$  ion



- PRESS THE MODE BUTTON 1 TIME TO TURN ON THE  $\oplus$  ION.
- PRESS THE MODE BUTTON ONE MORE TIME TO TURN ON THE  $\ominus$  ION.
- IF YOU PRESS THE ION ONE MORE TIME WHILE IT'S ON, THE DEVICE TURNS OFF.


\* PRESS AND HOLD FOR SOME SECONDS WHILE IN USE TO GENERATE VIBRATION.



## How to Use

5. Always hold the battery cover by hand, then brush it close to the scalp and brush it slowly as if you were brushing your hair.

\* Depending on the scalp product used, the direction or method of combing may vary.

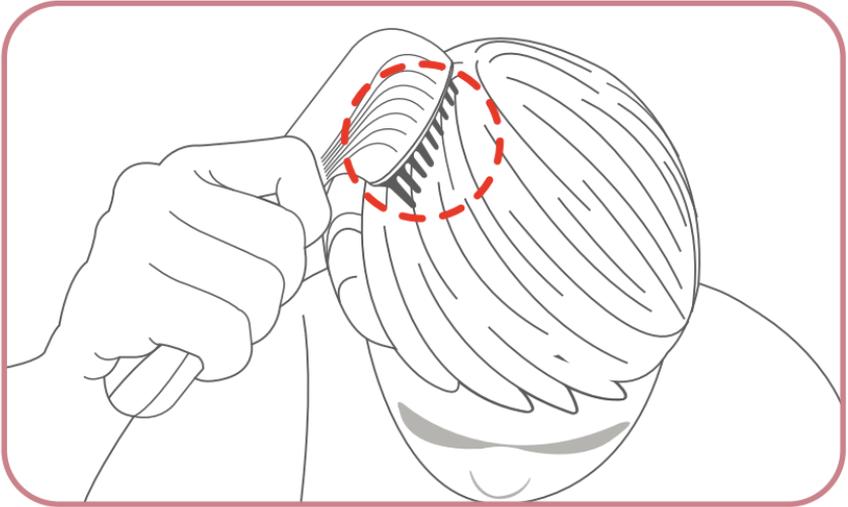

6. After use, always wash the brush with wet tissue, etc. and store it.


\* Make sure to wash it for the clean use next time.




# HOW TO USE

## Brush Exchange Method

1. Push the release button of the brush with your hand as shown to remove the brush from the device.

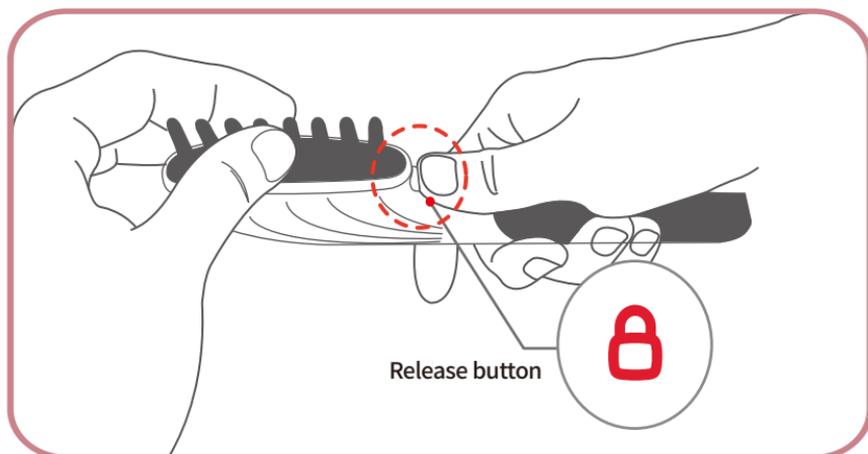

2. When you install a new brush, Hang the front first and press the back as shown in the picture to assemble the brush into the device.

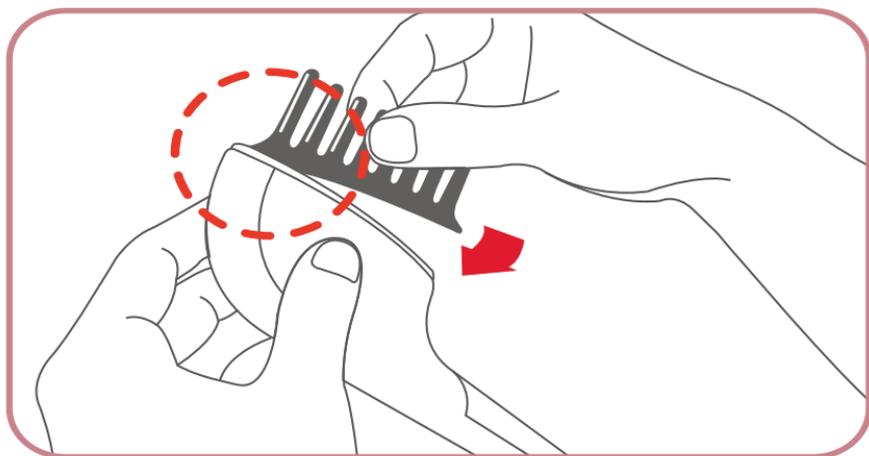

\* Brushes are consumables

## Cleaning and Storage Method

1. Make sure to wipe the brush clean with a wet tissue after using it. Wipe with a rough towel, etc. and the plating may come off.
2. When not in use for a long time, always keep the battery disconnected from the device.  
The fault caused by the leakage of batteries is consumer fault.
3. Keep it in a place that is well ventilated, moist, and free of harmful chemicals.
4. Please store it in a place without direct sunlight.
5. Consumer recommended storage temperature is 1~40°C and consumer recommended storage humidity is 20~60%.

## Q&A

- Q1. Pressing the Mode button does not turn on power.  
 → **The battery may have reached its end of life. Please replace it with a new battery.**  
 → **Use alkaline batteries for smooth use of the device.**
- Q2. I felt a prick in my scalp while using the device.  
 → **Please stop using it right away and contact us.**
- Q3. I can feel the vibration, but I don't know if the galvanic works well.  
 → **The battery cover must be held by hand to produce galvanic ions.**

## Product specification

- Model name: Reevaa H
- Manufacturing number: separate marking
- Date of manufacture: separate marking
- Weight (body): 65g
- Power: 3V (two AA batteries)
- Packaging Unit: 1 Set
- Date of preparation: [Revision] July 1, 2020  
(This product may be changed without notice to improve product performance.)

# A/S WARRANTY

- Model name: Reevaa H
  - Warranty period: One year from the day of purchase
  - Day of purchase:  
(Make sure to write down the day of purchase.)
  - Warranty might be void if the day of purchase is not written.)
- Purchased at:

When malfunctioning occurs within the warranty period while normally using the product in accordance with the user manual and cautions, free A/S will be based on the regulations for warranty.  
provided with charge in case the malfunctioning falls under the charged services in the regulations for warranty.

## **The followings will be charged even within the warranty period.**

- In case malfunctioning is caused by the mishandling of user (drop, water damage, breakage, etc)
- In case of exterior change (scratch, discoloring, etc) during user.
- In case malfunctioning is caused by not following the user manual or cautions.
- In case malfunctioning is caused by the unauthorized repair, disassembly, remodeling etc of the user, or by the carelessness of the user.
- In case malfunctioning or damage is caused by using separate products not provided by the seller, or by dropping the product while moving.
- In case this letter of warranty or the day of purchase is missing.
- In case malfunctioning occurs on the charger or expendables.
- In case malfunctioning is caused by natural disasters such as fires, earthquakes, flooding, thunderbolts, etc.